PERIO TRAY® THERAPY BY PERIO PROTECT

HOMECARE GUIDE

TERAPIA PERIO TRAY® DE PERIO PROTECT

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PerioProtect 

perioprotect.com

HG.2025/02

View our Patient
Resource Center

View our IFU video
on YouTube

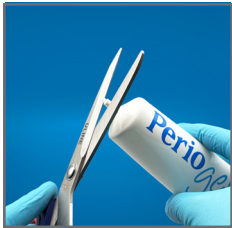
Figure A

Figure B

Figure C

INSTRUCTIONS FOR USE

1. Cut the tip of your Perio Gel® tube [A].
2. Place or hold the Perio Tray® open side up [B].
3. Apply a small, thin line across the tooth indentations [C].
4. Place the tray into your mouth so the tray covers all the teeth and fully seats in place. Spit out any excess medication.
5. Wear your trays for the time prescribed by your doctor. (Typically, 15 minutes per session)
6. After use, remove the trays. Brush your teeth to remove any excess medication.
7. Use your toothbrush to clean the trays and rinse them thoroughly with cold water. **DO NOT USE HOT WATER.**
8. Gently towel-dry your trays and return them to the storage case.

It is very important that you continue to brush and floss your teeth. If you experience any adverse reactions during Perio Tray® Therapy, discontinue use and contact your dentist.

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PERIO TRAY® WEARING INSTRUCTIONS

FOR PATIENTS DIAGNOSED WITH PERIODONTITIS

- 3 times a day for 15 minutes
- 2 times a day for 15 minutes
- 1 time a day for 15 minutes

FOR PATIENTS DIAGNOSED WITH GINGIVITIS

- 2 times a day for 15 minutes
- 1 time a day for 15 minutes

FOR MAINTENANCE

- 2 times a day for 15 minutes
- 1 time a day for 15 minutes

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